


## **OTF 5000**

Microtome Cryostat

## The Anatomy of the OTF5000

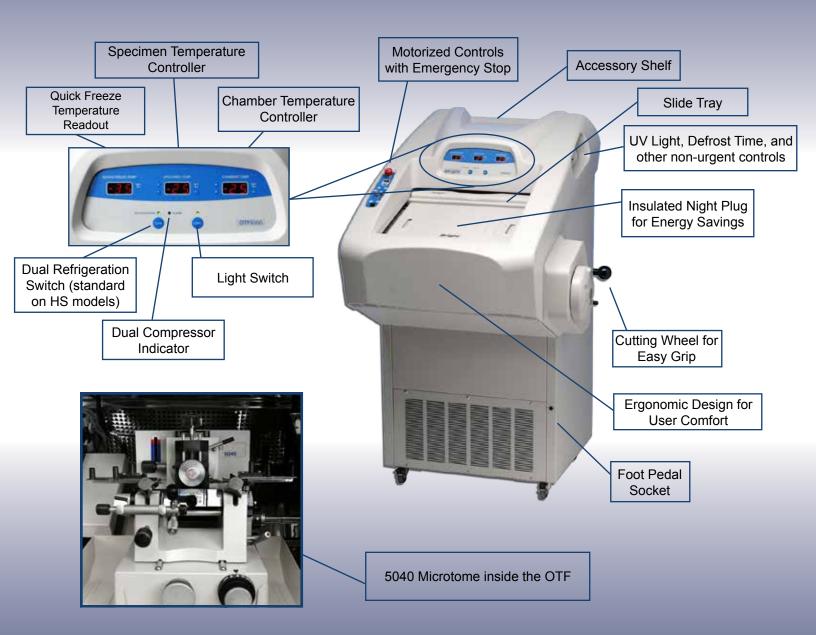

www.hackerinstruments.com www.tidylabstore.com 1-800-4-HACKER HACKERLAB@AOL.COM



### OTF 5000 Packages

OTF/LS-002	-35° operation with heavy duty single compressor, auto-defrost and buffered quartz tir freezer digital temperature displays; Model 5040 rotary-retracting microtome with 1/60µ and standard knife block.  Specimen Temperature control:			
Especially geared for brain sectioning and other difficult-to-section tissue such as breast.	<ul> <li>Allows rapid changes to be made to the specimen temperature, above the chamber temperated.</li> <li>Creates a positive temperature difference between the knife and the specimen so the will always be delivered onto a surface colder than itself.</li> <li>Ideal when different tissues need to be cut at different temperatures.  Object orientation with fine adjustment mechanism:  Allows for accurate orientation of the specimen in both X &amp; Y planes.  Includes both macro and micro adjustment knobs.  Firm clamping system for error free sections.  Motorized advance and rewind for rapid control of head position:  Two-speed advance/rewind provides slow, rapid, or incremental advances.  Auto reset returns cutting head to start position quickly.  Also Features: Adjustable anti-roll plate for cutting difficult specimen such as brain.</li> </ul>			
	Cat. No. OTF 5000/LS-002 \$3	32,500.00		
OTF/LS-004	Includes all features of the LS-002 along with:  Microtome with Quick Release Feather Blade Holder in place of Standard Knife Block. All other features as listed above.			
	Cat. No. OTF 5000/LS-004 \$3	32,950.00		
OTF/HS-002  Suitable for any tissue or material sectioning  For advanced research and extensive experimentation with automated sectioning	Cat. No. OTF 5000/LS-004 \$3  Includes all features of the LS Models, along with:	just ng velocities sue sizes		
Suitable for any tissue or material sectioning  For advanced research and extensive experimentation with automated	<ul> <li>Includes all features of the LS Models, along with:         <ul> <li>Dual Refrigeration System:</li> </ul> </li> <li>Individually selectable dual compressors for 24/7/365 reliability</li> <li>Ability to reach -45° C chamber temperature with both compressors</li> <li>Advanced sharing algorithm to extend individual compressor life; reliable logic to adj compressor cycles to maintain selected temperature even upon compressor failure Robomatic Motor Drive for automated or manual sectioning:         <ul> <li>Provides hands-free operation for ergonomic use as well as extremely smooth cuttin for continuous or single sections</li> <li>Advance cutting window to adjust speeds at different cutting locations for various tiss</li> <li>Variable speed capability to match varying tissue textures</li></ul></li></ul>	just ng velocities sue sizes		
Suitable for any tissue or material sectioning  For advanced research and extensive experimentation with automated	<ul> <li>Includes all features of the LS Models, along with:         <ul> <li>Dual Refrigeration System:</li> </ul> </li> <li>Individually selectable dual compressors for 24/7/365 reliability</li> <li>Ability to reach -45° C chamber temperature with both compressors</li> <li>Advanced sharing algorithm to extend individual compressor life; reliable logic to adj compressor cycles to maintain selected temperature even upon compressor failure Robomatic Motor Drive for automated or manual sectioning:         <ul> <li>Provides hands-free operation for ergonomic use as well as extremely smooth cuttin for continuous or single sections</li> <li>Advance cutting window to adjust speeds at different cutting locations for various tiss</li> <li>Variable speed capability to match varying tissue textures</li></ul></li></ul>	just ng velocities sue sizes nter and		

Standard Accessories: 3 plain object holders, 18cm "C" profile knife or pack of disposable blades, knife cleaning brush, 3 Allen keys, spare fuses, Cryo-M-Bed, spare anti-roll plate, debris tray, utility shelves, debris plug, low temperature oil, operating manual

Instrument Dimensions:
49" (H) x 33" (D) x 30-35" (W) (manual/motor drive);
Weight 265 lb min / 420 lb max

# OTF 5000 Ancillaries & Replacements

	Cat. No.	Description	PRICE USD
	10 70 12	0.2/12 μm	\$1,812.00
	10 70 13	0.5/30 μm	\$1,590.00
Inter-changable	10 70 14	1/60 μm	\$2,053.00
advance	10 70 50	2/120 μm	\$2,168.00
movements	10 70 55	3/180 μm	\$2,205.00
	10 70 56	5/300 μm	\$3,415.00
	10 70 07	5 mm diameter	\$107.00
Plain object	10 70 08	22 mm diameter	\$86.00
holders; use directly	10 70 09	37 mm diameter	\$202.00
in microtome or "O" ancillary (10 70 18)	10 70 10	50 mm diameter	\$252.00
andmary (10 70 10)	10 70 11	50 x 50 mm diameter	\$483.00
Plain object holders,	10 70 08-1	22 mm longer length	\$86.00
longer length	10 70 09-1	37 mm longer length	\$202.00
	10 70 201	Universal holder for 5, 23, 37, & 50 mm object holders	\$1,340.00
	10 70 21	5 mm diameter, aluminum	\$140.00
Orientable object holders	10 70 22	22 mm diameter, stainless steel	\$89.00
noiders	10 70 23	37 mm diameter, stainless steel	\$180.00
	10 70 24	50 mm diameter, stainless steel	\$295.00
		Order 5 or more of one size orientable object	
		holder for 10% discount	
Tubular object holders. Not for use with	10 70 26	22 mm diameter	\$140.00
"O" ancillary (10 70 18)	10 70 30	Object holder - 50 x 70 mm	\$445.00
For large sectioning,	10 70 32	Modification to standard anti-roll device for 70 mm sections	\$2,525.00
e.g. brain, whole	10 70 34	Knife extension plate to hold 50 x 70 mm sections	\$556.00
small animals	10 70 38	Modification to microtome - extension of knife block excursion	\$315.00
	10 70 45	Trough and knife conductor, right side	\$390.00
	10 70 47	Trough and knife conductor, left side	\$390.00
Microtome	13 27 18	H/I "Cryo-C" Cryostat knife, Profile "C", 22º, 18 cm	\$540.00
Knives	13 20 16	Profile C, 16 cm	\$450.00
	13 20 18	Profile C, 18 cm	\$510.00
For sectioning	13 00 80	Tungsten Carbide Edge, Profile "D" - 16 cm	\$2,850.00
calcified tissue	13 00 82	Tungsten Carbide Edge: Profile "C" 22º - 22/16 cm	\$3,200.00
For aloes knives	10 10 05	Glass knife holder for "Ralph" knives	\$1,750.00
For glass knives	10 70 15	Pyramid Glass Knife Holder, (6-12 mm glass) fully orientable	\$2,600.00
	10 75 40	SuperSlip™ Glass Coverslips 24 x 40 (10oz)*	\$115.00
	10 75 50	SuperSlip™ Glass Coverslips 24 x 50 (10oz)*	\$115.00
Coverglass	10 75 60	SuperSlip™ Glass Coverslips 24 x 60 (10oz)*	\$115.00
		*Special pricing considerations for coverglass are available, contingent on volume	
		are avaliable, continuent on volume	

All prices subject to change without notice


## OTF 5000 Ancillaries & Replacements

	Cat. No.	Description	PRICE USD
"Magnacut" special disposable knife system	10 70 63	"Magnacut" magnetic knife holder	\$1,750.00
	11 70 64	Magnetic knife holder with Magnaplate anti-roll	\$1,975.00
	11 70 50	Magnetic knife holder with standard anti-roll	\$2,100.00
	10 70 62	Magnaplate anti-roll guide plate (part of 10 70 64) for use with disposable knives	\$105.00
	11 55 10	Disposable Magnacut blades, packet of 10	\$150.00
	11 70 53	High/low profile blade holder (uses 11 70 52 or 11 70 58)	\$1,550.00
	11 70 56	High Profile blade holder (uses 11 70 58 only)	\$1,550.00
	11 70 54	Low Profile blade holder (uses 11 70 57 only)	\$1,550.00
	11 70 55	Heavy duty blade holder w/ extended handle (uses 11 70 57 only)	\$1,550.00
Miscellaneous	11 70 58	H/l Disposable high profile blades, pkg of 60 ea	\$183.00
	11 70 52	H/I Disposable low profile blades, pkg of 50 ea	\$180.00
	11 70 57	H/I Disposable heavy duty blades, pkg of 35 ea	\$204.00
	11 70 65	Lab Tek object holder, adaptor, for quick-freezer, set of 2	\$410.00
	11 70 82	Set of 2 shelves	\$570.00
For storage of	11 70 18	Quick release object holder for microtome clamp, complete with (50) 22 mm object discs	\$2,425.00
samples of LN2	11 70 19	Extra set of (50) discs	\$595.00
	10 70 72	Shaft, 72 mm long, for anti-roll plates	\$42.00
	10 70 69	Anti-roll plate, 37 mm long	\$170.00
Anti-Roll Plates	10 70 66	Anti-roll plate, 50 mm long	\$170.00
	10 70 67	Anti-roll plate, 76 mm long	\$170.00
	10 70 68	Anti-roll plate, 93 mm long	\$170.00
	11 70 74	Fluorescent tube	\$26.00
	11 70 76	Cyro-M-Bed, 113 ml (4oz)	\$57.00
Replacements	11 70 25	Low temperature oil, amber	\$78.00
	11 40 49	Installation & Operation Manual	\$95.00
	10 70 19	Service Manual	\$105.00

All prices subject to change without notice

Cryostat ancillaries, replacement parts and consumables available at www.tidylabstore.com

#### OTF 5000 Knives & Holders

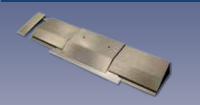

#### Lever Release Disposable Blade Holder

The disposable blade holder fits into the existing knife block holder for use of disposable and feather blades. It is very easy to install and the entire length of the blade can be utilized. Ideal for all soft tissues, and very good for brain tissue



#### Standard Knife Block

The standard knife block holder permits the use of both C and D profile steel knives, as well as tungsten carbide knives. The standard knife block provides easy adjustment of the knife angle and can easily be changed with the quick release feather blade holder



#### Quick Release Feather Blade Holder

The Quick Release Feather Blade holder allows the use of the entire blade length and allows the angle of the blade to be adjusted without adjusting the blade position. Feather blades are ideal for all soft tissue, and a very good solution for brain sectioning



#### Magnacut Disposable Blade Holder

Long-lasting Magnacut blades provide the flexibility to use disposable blades and increases the time between replacements. The blades are held by a strong magnet providing security and precision. The blade adjustment tool permits the entire length of the blade to be utilized. Magnacut blades can be used for all soft tissues except brain - they can also be used on tougher tissues, but not bone.



#### Hacker "Cryo-C" Microtome Knife

Specifically designed for the OTF Cryostat. 18cm "C" profile knife with a 22° cutting angle



#### Hacker Disposable Knives

Available in high and low profile, heavy duty and "Select" grade



## www.hackerinstruments.com www.tidylabstore.com

1132 Kincaid Bridge Road • PO BOX 1176 Winnsboro, SC 29180 • (803) 712-6100



# A Perfect Match Bright OTF 5000 Microtome Cryostat *from* Hacker

# The **Hacker** Family of Cryostats



2212 Micro
Cryostat
Portable BenchTop Micro-Crypstat



8250 & 9400
Electro-Linear Sledge Microtome
Full Body Cryostats



Clinicut
Free Standing
Mobile Clinical
Microtome Cryostat

Visit us online at www.hackerinstruments.com for more information