

H/I Hacker CellSpin CytoCentrifuge

a multi- functional system for producing mono layer cell preparations

Advantages & Benefits

- ◆ The first centrifuge that grows with your laboratory's requirements
- ◆ An affordable 4 to 12-position centrifuge
- ◆ With the "Double Cellfunnel" you can double the number of preparations.
- ◆ You work with the well-known monolayer filter card system
- ◆ Acquire outstanding results routinely for thin-layer cytology

Range of applications

Suitable for cytology, urology, microbiology, hematology, immunohistochemistry, virology, oncology - for small and large volumes!

The cells are deposited on an area of 6 mm or 22 x 15 mm.



Safety features

Metal case and lid • safety catch to prevent lid from falling • lid lock emergency release* • motor overheating protection • unbalance cut-off • simple rotor change • automatic rotor recognition

Specifications

- Programmable up to 99 minutes
- Up to 3 programs can be saved
- The speed is selectable between RPM or RCF from 100 to 6,000 rpm or 100 to 18.000 rpm.
- The acceleration can be set to fast, medium or even slow in order to protect the cells
- The stopping times are programmable for gentle cell treatment
- Alarm function indicates unbalance

The CellSpin comes in five versions:

CellclipRotor	CellSpin I-4	CellSpin I-8	CellSpin I-12	CellSpin II-12	CellSpin II-12 R > chilled
	4 positions	8 positions	12 positions	12 positions	12 positions
Clips	4 cell clips	8 cell clips	12 cell clips	12 cell clips	12 cell clips
Order No.	JC 100230/4	JC 100230/8	JC 100230/12	JC 200230/12	JC 200230/12R
	Incl. starter kit for 200 tests <input type="checkbox"/> Double the volume with the double Cellfunnel!				

Reusable and disposable sample chambers



JC320 Cellfunnel disposable sample chamber, 100/500 units
JC304 Cellfunnel reusable chamber, 12 units



JC305 Filter cards for Cellfunnel SC304, 200 units
JC310 CytoSlide - single cycle, uncoated, 100 units
JC311 Cyto Slide - single cycle, coated, 100 units



JC306 Double Cellfunnel reusable chamber, 12 units
JC307 Filter cards for double Cellfunnel JC306, 500 units
JC321 Double Cellfunnel disposable chambers, 100/500 units
JC309 Cyto Slide - double cycle, uncoated, 100 units
JC312 Cyto Slide - double cycle, coated, 100 units




JC313 ECO funnel reusable chamber, 12 units
JC315 ECO seal for ECO funnel, 100 units
JC322 ECO funnel disposable chamber, 100/500, units
JC314 ECO micro slide carrier, coated, 100 units

Cellcliprotor



JC101 The optionally obtainable, sealed rotor handles up to 24 samples per run. It protects the operating staff when loading, using and unloading into the safety chamber. The rotor is autoclavable.

Cellcheck

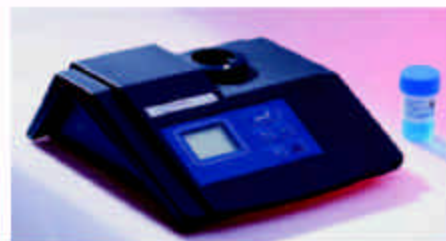

JC339 You can determine the optimal number of cells for producing preparations using the Cell check.

JC390/T CellSpin replacement cell clip
JC390/S Shandon compatible cell clip

HACKER Instruments & Industries Inc

1132 Kincaid Bridge Road • P.O. Box 1176
Winnsboro, South Carolina 29180, USA
Tel: (803) 712-6100 • Fax: (803) 712-6116
e-mail: hackerlab@aol.com
www.hackerinstruments.com